
CHŁODZIARKA FARMACEUTYCZNA PHARMAZEUTISCHE KÜHLSCHRÄNKE PHARMACEUTICAL REFRIGERATORS


850066, 850350

**PL**
3-6**DE**
7-10**EN**
11-14

85066

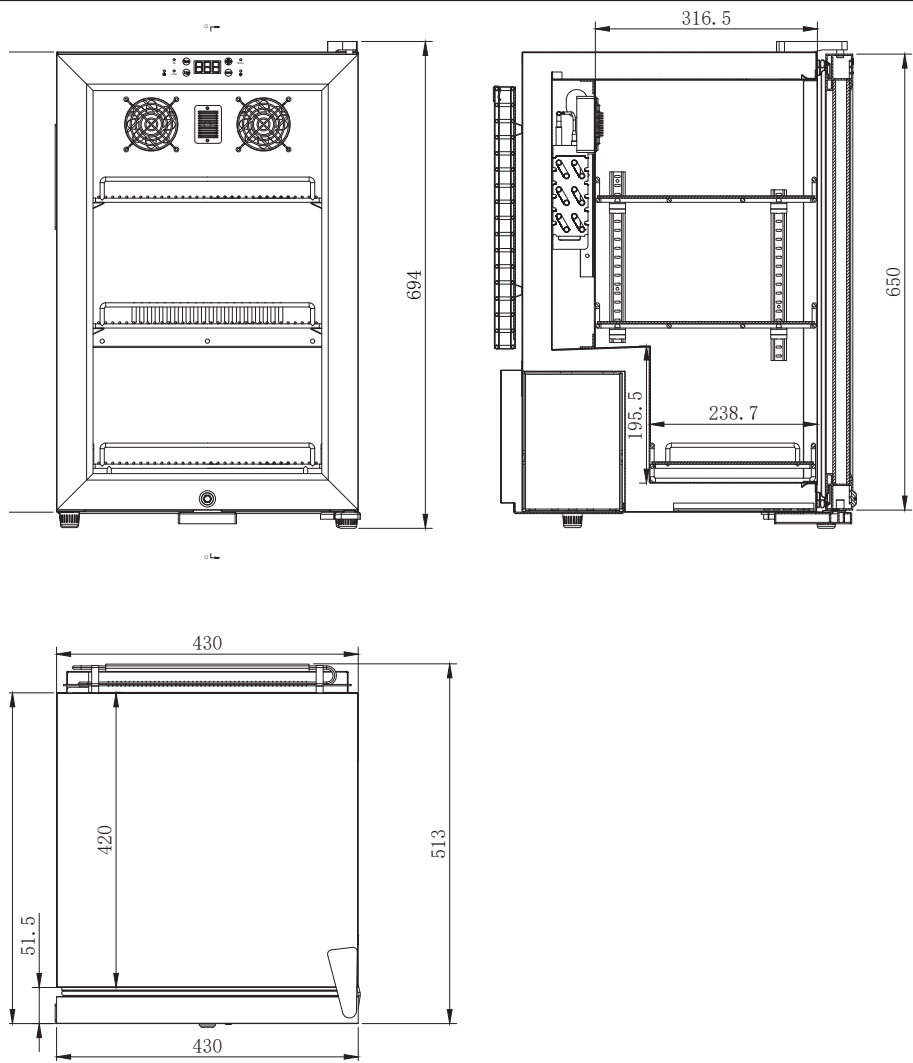850350

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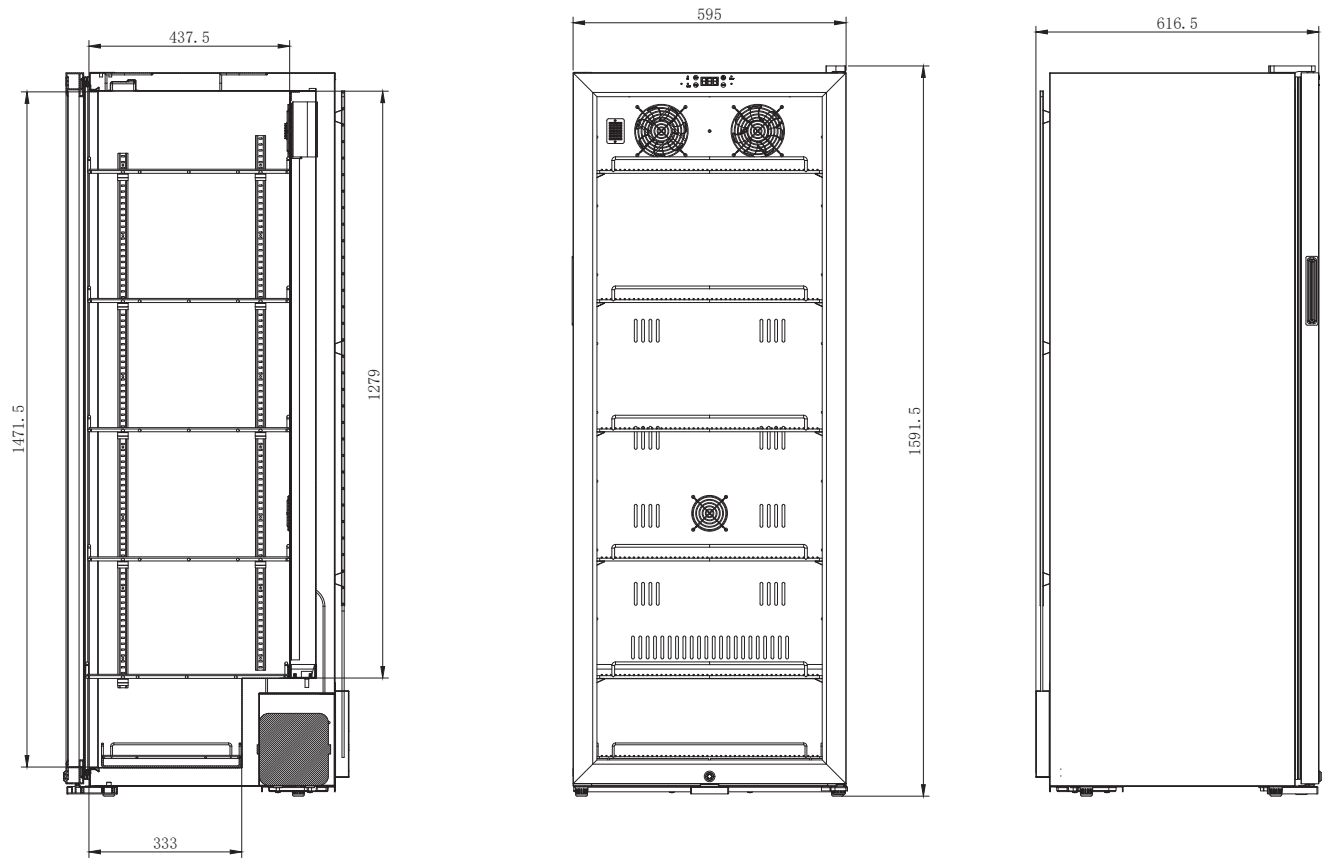Thank you for the purchase of our product. Please read this manual carefully before starting to operate the appliance. All rights reserved. No part of this manual may be reproduced by photocopy, printing or in any other way without prior consent from the manufacturer.

Photographs and drawings are furnished by way of example and can vary from the purchased appliance.

NOTE: Keep the manual in a safe place accessible for the personnel.

The manufacturer reserves the right to change technical parameters of the appliance at any time without advance notice.

1. SAFETY INSTRUCTIONS

- Incorrect operation and improper use can result in serious damage to the appliance or personal injury.
- Use this appliance solely for its intended use as described in this manual.
- The manufacturer cannot be held liable for any damage caused by incorrect operation and improper use of the appliance.
- During operating keep the appliance and the plug away from water or other liquids. In case the appliance, as a result of carelessness, has been immersed in water or flooded by water immediately unplug the appliance from the wall socket and prior to further use contact a qualified technician for examination.
- Failure to observe instructions enclosed in this manual may result in the hazard of injury or death.
- Do not let unauthorized opening of the housing of the appliance.
- Do not insert any foreign objects into the housing of the appliance.
- Do not touch the plug and the power cord with damp hands.
- Regularly check the condition of the plug and the cord. If the power cord or the plug is damaged they should be replaced by the Authorized Technical Service personnel.
- In case the appliance is dropped or damaged in any other way before turning it on always contact the Authorized Technical Service for examination and repair.
- Any repairs and work on the appliance should be carried out solely by Authorized Technical Service, never by the end user as unauthorized tampering could prove highly dangerous for the user.
- Do not let the cord hang over or touch hot surfaces or sharp edges. Keep the cord away from naked fire. When disconnecting the appliance from the wall socket, pull on the plug, never on the cord.
- Do not let the cord (or extension cord) be inadvertently pulled on or tripped over.
- Check the operation of the appliance during usage.
- This appliance can be used by children in the age of eight and more or by persons with reduced physical, sensory or mental capabilities or lack of knowledge and experience in operating of the appliance, if they have been given supervision or instruction concerning use of the appliance in a safe way by a person responsible for their safety and they understand the hazards involved. Children should not play with the appliance. Children should not perform cleaning or maintenance of the appliance without supervision.
- Always switch off the appliance and unplug the power cord from the wall socket when the appliance is not in use and prior to cleaning operations. The disconnected plug should be sited in a visible place.
- **Warning:** The appliance's power will remain ON unless unplugged.
- Turn off the appliance before pulling the plug out of the wall socket.
- Never pull the appliance on the cord.
- It is prohibited to store any explosive substances such as aerosol / flammable gas cans in the appliance.
- **WARNING:** Do not block or cover any ventilation openings in the equipment casing or in the built-in structure.
- **WARNING:** Do not use mechanical or any other means to accelerate the defrosting process other than recommended by the manufacturer.
- **WARNING:** Protect the refrigeration system from damage.
- **WARNING:** Do not use any type of electrical equipment other than recommended by the manufacturer inside the food storage chamber.

2. DESCRIPTION AND INTENDED USE OF THE APPLIANCE

- A pharmaceutical refrigerator is a specialized device designed to store medicinal products, vaccines, reagents, biological samples, and other temperature-sensitive preparations under strictly controlled conditions.
- forced air circulation
- temperature controller
- automatic defrosting

3. TECHNICAL SPECIFICATION

Model	850066	850350
Power supply (V/Hz)	230 / 50	
Current (A)	0,4	0,35
Power (W)	135	250
Temperature range (°C)	2 / 8	2 / 8
Dimensions (WxDxH) (mm)	strona 2	
Climate class *	4	4
Weight (kg)	35	82
Refrigerant / quantity (g)	R600a / 23	R600a / 50
Capacity (l)	63	300

* 4 - to +30°C / 55% relative humidity

The device contains R290 or R600a coolant - natural, ecological, flammable gas.

Take special care when transporting and installing the device to avoid damaging the cooling circuit.

If a leak is revealed, refrain from using devices/items that could cause sparks or fire and ventilate the room. To avoid the formation of a flammable gas mixture in the event of a leak from the refrigeration system, install the device in a room with area of 1m³ per every 8 g of the coolant. The amount of coolant is given in „TECHNICAL DATA“ and on the nameplate.

4. TRANSPORT AND STORAGE

The appliance is protected with a protective film, placed on a wooden pallet, packed in a cardboard box and secured with stripes.

- packaged appliance should be stored in a covered warehouse at an ambient temperature of 0 °C / + 55 °C and humidity of 30 - 95 %;
- it is forbidden to stack appliances;
- read the data on the appliance rating plate. In case of damage or loss of the rating plate it should be replaced immediately;
- do not unscrew the safety guards;
- after unpacking the appliance, carefully remove the protective film and other elements protecting the appliance during transport;
- while lifting the appliance with a lift, position the lift on the outside of the frame;
- transport the appliance and store it exclusively in a vertical position. Never place the appliance on its side or upside down, because it can damage the cooling system (not subject to warranty repair);
- after transporting or moving the appliance wait 2 hours before starting it;
- do not place the appliance in the vicinity of flammable products, outdoors or in a place subject to direct sunlight exposure.

5. INSTALLATION

5.1. Positioning of the appliance

The appliance should be placed in a dry, adequately ventilated room. In order to ensure proper work do not place the appliance near a heat source and do not expose it to direct sunlight. The optimal temperature range in the room in which the appliance operates is between: + 16 °C and + 35 °C.

Avoid placing the appliance in rooms with a high chlorine and acids content in the air (e.g. at a pool), as this may cause corrosion of the surface made of stainless steel.

The cooling cabinet is prepared to be positioned at a wall as a free standing appliance.

Important: The cooling cabinet should be installed in a place with proper ventilation ensuring free air circulation around the appliance. The correct distance from the rear panel of the appliance to a wall is ensured by special elements placed at the back of the cabinet. Do not place the appliance in the vicinity of flammable products, in the open air or in a place exposed to direct sunlight.

Level the device using the adjustable feet.


5.2. Electrical connection

The appliance should be connected to an AC power source.

Before connecting the appliance to the mains, make sure that the voltage and frequency values given on the rating plate of the appliance correspond with the voltage and frequency values at the place of installation.


The electrical outlet should be properly earthed and equipped with a differential current protection device in accordance with legal provisions in force.

6. CONTROL PANEL

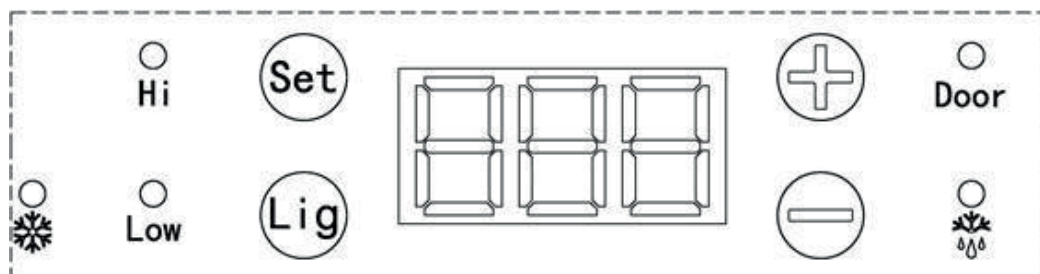6.1. Button functions

Button / indicator light	Function
	Turn the unit on/off.
	Turn the interior lighting of the refrigeration chamber on/off.

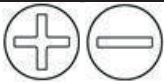Button / indicator light	Function
	Increase/decrease the temperature.
	Cooling indicator.
	High/low temperature alarm indicator inside the refrigeration chamber.
	Defrost indicator.
	Door open alarm indicator.

6.2. Button functions

1. Setting user parameters

Press and hold the button  for 5 seconds, then use the  or  buttons to adjust the set temperature. After about 5 seconds, the chamber temperature will be displayed. Each press of the temperature adjustment buttons increases the temperature by 0.1 °C.

If no button is pressed within 30 seconds in any setting mode, the newly set temperature will be saved automatically.

Light switch  – each time you press it, the light turns on; pressing it again turns it off. After 10 minutes of inactivity, the lighting turns off automatically because certain biological products, chemical samples, pharmaceutical preparations, cosmetics, etc., are light-sensitive, and prolonged exposure to light can affect their properties. Therefore, the built-in lighting emits a cool, soft light that provides illumination without negatively affecting the chemical properties of the stored products.

2. Defrosting

The appliance is equipped with an automatic defrosting system. When the ", the device is currently defrosting.

7. TROUBLESHOOTING

Every time the power is turned on, all components operate, all indicator lights come on, and the system performs a 10-second self-test. The compressor starts up 3 minutes after the self-test is successfully completed.

7.1. Alarm code, root cause analysis, and corrective actions

Alarm code	Reason	Solution
The "Hi" or "Low" indicators flash alternately with the temperature display. An audible alarm sounds.	The temperature inside the chamber is higher/lower than the set value, and the refrigerator door remains open for an extended period.	Close the door and check the set temperature. The alarm will turn off automatically.
"Door" indicator: Flashes alternately with the temperature display. An audible alarm sounds.	The door is not closed or is not closed properly.	Close the door, and the alarm will turn off automatically.
'E1' or "E2" flash alternately with the temperature display and an audible alarm.	Failure of the cooling system temperature sensor.	Contact an authorized service center or dealer.
"E3" or "E4" - alternates with the temperature display. An audible signal sounds.	Failure of the defrost temperature sensor.	Contact an authorized service center or dealer.
"E5" – alternates with the temperature display. An audible signal sounds.	Failure of the defrost system.	Contact an authorized service center or dealer.
"E6" - alternates with the temperature display. An audible signal sounds.	Failure of the compensation heater. Communication between the power supply board and the main board is disrupted.	Contact an authorized service center or dealer.

IMPORTANT! If any of the above errors appear on the display, contact customer service immediately and provide the error code.

8. CLEANING AND MAINTENANCE

Note: Before starting any maintenance procedure, always disconnect the appliance from the power supply.

The appliance should be washed with water and dishwashing liquid. It is not recommended to clean the appliance with detergents, because they can damage the walls of the cabinet. Plastic parts of the appliance can be washed with water at the maximum temperature of 85 °C.

During daily cleaning of the appliance, remember to thoroughly clean the gasket from crumbs and sticky substances, because they can damage it.

The ventilation grate of the compressor (located on the back of the appliance) should be cleaned regularly, as dirt and dust may block free air circulation.

9. PROCEDURE IN THE EVENT OF POWER FAILURE


In case the appliance is disconnected from the power supply (both due to a failure or due to accidental removing of the plug from the electric socket), the element controlling the work of the compressor will cause that the restarted appliance will operate according to the current temperature control knob setting.

10. DISPOSAL OF THE OLD APPLIANCE

Information for users on the rules concerning proper disposal of electrical and electronic equipment at the end of its life.

- **Leave your old equipment in the shop where you purchase a new appliance.** Each store is obliged to accept old equipment free of charge if you buy there the new equipment of the same type and in the same quantity. The only condition is to deliver the old equipment to the store at your expense.
- **Dispose of your old equipment at a collection point.** Information about the location of the nearest collection point can be found on the local website or the notice board of the municipal office.
- **Leave the equipment at a service point.** If the repair of the equipment is unprofitable or impossible for technical reasons, the service is obliged to accept this appliance free of charge.
- **Return the used equipment at the end of its life without leaving your home.** If you do not have the time or ability to deliver your equipment to the collection point, you can use the services of specialized companies.
- **Remember! Do not dispose of used equipment together with regular waste.**

An illegal or incorrect disposal of the product will lead to high penalties being imposed.

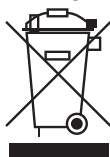The symbol of a crossed-out rubbish bin on the product, its packaging or its instruction manual means that at the end of its life the product must not be disposed of with regular waste. The user is obliged to deliver their used equipment to a designated collection point for proper processing.

11. GUARANTEE

The dealer of the appliance is liable under the guarantee.

The guarantee does not cover any damage resulting from the build-up of limescale inside the appliance. The guarantee also does not cover: damage caused by external forces such as atmospheric discharge, voltage change, incorrect setting of electrical voltage, power supply from unsuitable power outlet, any mechanical, thermal, chemical damage and defects resulting from them.

Furthermore, the guarantee does not cover the replacement of exhaustible parts, namely: light bulbs, rubber elements, heating elements damaged by limescale, glass and other elements that are exposed to any form of wear or usage, for example; burners, rubber units as well as any breakdowns or mechanical damage of the appliance elements caused by incorrect storage, transport, maintenance, installation or operation.

Stalgast Sp. z o.o.
ul. Ostrobramska 75C, lokal 6.02, 04-175 Warszawa
tel.: 22 517 15 75 fax: 22 517 15 77
www.stalgast.com email: stalgaſt@stalgaſt.com

• DE •
STALGAST GmbH
Mary-Somerville-Str. 6,
28359 Bremen;
Tel.: +49 421 9898066-1
stalgaſt@stalgaſt.de
www.stalgast.de

• EN •
Tel.: +48 22 509 30 77
export@stalgaſt.com
www.stalgast.eu

• FR • ES •
Tel.: +48 22 509 30 55
export@stalgaſt.com
www.stalgast.eu